



OZP 558



OZP 555



Trinocular head

## LAB LINE

Professional and powerful – thanks to its extremely large magnification range, strong illumination and first-class optics

### Features

- The KERN OZP stereo zoom microscope stands out through its above-average magnification range and its robust shape which is also ergonomic, it enables effortless, simple working over a period of several hours.
- The KERN OZP series is available as a strong, continuously adjustable 3 W LED reflected and transmitted light variant for the very best illumination of your sample or as a variant without illumination.
- With its large working distance, an extra large field of view and brilliant resolution, the KERN OZP provides sharp, high-contrast and colour-true images.
- The extremely large, continuously adjustable magnification range from 6 to 55 times magnification means that you can work quickly and effectively.
- There is a choice of two binocular models as well as two trinocular models for connecting a camera for documentation purposes and for quality reports.
- The pillar stand is particularly flexible due to its variable and sturdy adjustment mechanism and therefore enables ergonomic working procedures.
- A large selection of eyepieces, (universal) stands, a darkfield kit, external illumination units as well as auxiliary objectives and more are available as accessories.
- A protective dust cover, eye cups, as well as multi-lingual user instructions are included in the scope of delivery.
- A C-mount adapter is required to connect a camera to the trinocular version. You can select this adapter from the following model outfit list.
- Please find detailed information in the following charts.

### Scope of application

- Zoology and botany, quality control, electronics and semiconductor industry, assembly and repair

### Applications/Samples

- Samples with focus on three-dimensional impression, zoom with variable magnification (depth, thickness), e.g. insects, seeds, circuit boards, components

### Technical data

- Optical system: Greenough optics
- Brightness adjustable (separate)
- Tube: 35° inclined
- Magnification ratio: 9,2:1
- Interpupillary distance: 52 – 76 mm
- Diopter adjustment: Both-sided
- Overall dimensions W×D×H  
330×285×470 mm
- Net weight approx. 4,5 kg

#### STANDARD



#### OPTION


























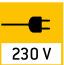








Model	Standard configuration						
	Tube	Eyepiece	Field of view mm	Objective Zoom	Stand	Illumination	
<b>KERN</b>							
<b>OZP 555</b>	Binocular	HSWF 10×/ø 23 mm	ø 38,3 – 4,2	0,6× – 5,5×	Pillar style	–	
<b>OZP 556</b>	Binocular	HSWF 10×/ø 23 mm	ø 38,3 – 4,2	0,6× – 5,5×	Pillar style	3 W LED (incident); 3 W LED (transmitted)	
<b>OZP 557</b>	Trinocular	HSWF 10×/ø 23 mm	ø 38,3 – 4,2	0,6× – 5,5×	Pillar style	–	
<b>OZP 558</b>	Trinocular	HSWF 10×/ø 23 mm	ø 38,3 – 4,2	0,6× – 5,5×	Pillar style	3 W LED (incident); 3 W LED (transmitted)	

Eyepiece	Specifications – Objectives						
	Magnification	Standard 1,0×	Auxiliary objectives				
			0,37×	0,5×	0,7×	1,5×	2×
HSWF 10×	Total magnification	6× – 55×	2,22× – 20,4×	3× – 27,5×	4,2× – 38,5×	9× – 82,5×	12× – 110×
	Field of view mm	∅ 38,3 – 4,2	∅ 99,1 – 10,8	∅ 76,7 – 8,4	∅ 54,8 – 6	∅ 25,6 – 2,8	∅ 19,2 – 2,1
SWF 15×	Total magnification	9× – 82,5×	3,33× – 30,5×	4,5× – 41,25×	6,3× – 57,75×	13,5× – 123,75×	18× – 165×
	Field of view mm	∅ 28,3 – 3,1	∅ 76,6 – 8,35	∅ 56,7 – 6,2	∅ 40,5 – 4,4	∅ 18,9 – 2,1	∅ 14,2 – 1,5
SWF 20×	Total magnification	12× – 110×	4,44× – 40,7×	6× – 55×	8,4× – 77×	18× – 165×	24× – 220×
	Field of view mm	∅ 23,3 – 2,5	∅ 63,1 – 6,9	∅ 46,7 – 5,1	∅ 33,3 – 3,6	∅ 15,6 – 1,7	∅ 11,7 – 1,3
SWF 30×	Total magnification	18× – 165×	6,66× – 61,1×	9× – 82,5×	12,6× – 115,5×	27× – 247,5×	36× – 330×
	Field of view mm	∅ 15 – 1,6	∅ 40,5 – 4,4	∅ 30 – 3,3	∅ 21,4 – 2,3	∅ 10 – 1,1	∅ 7,5 – 0,8
Working distance		108 mm	275 mm	195 mm	145 mm	50 mm	35 mm

Model outfit		Model KERN				Order number	
		OZP 555	OZP 556	OZP 557	OZP 558		
Eyepieces (30,0 mm)	HSWF 10×/∅ 23 mm	✓✓	✓✓	✓✓	✓✓	OZB-A5503	
	SWF 15×/∅ 17 mm	○○	○○	○○	○○	OZB-A5504	
	SWF 20×/∅ 14 mm	○○	○○	○○	○○	OZB-A5505	
	SWF 30×/∅ 9 mm	○○	○○	○○	○○	OZB-A5506	
	HSWF 10×/∅ 23 mm (reticule 0,1 mm)	○	○	○	○	OZB-A5512	
	SWF 15×/∅ 17 mm (reticule 0,05 mm)	○	○	○	○	OZB-A5513	
	SWF 20×/∅ 14 mm (reticule 0,05 mm)	○	○	○	○	OZB-A5514	
Achromatic auxiliary objectives	0,37× only in combination with a universal stand	○	○	○	○	OZB-A5611	
	0,5×	○	○	○	○	OZB-A5612	
	0,7×	○	○	○	○	OZB-A5613	
	1,5×	○	○	○	○	OZB-A5615	
	2,0×	○	○	○	○	OZB-A5616	
	Soldering protection lens	○	○	○	○	OZB-A5614	
C-Mount	0,3× (focus adjustable)			○	○	OZB-A5701	
	0,5× (focus adjustable)			○	○	OZB-A5702	
	1,0× (focus adjustable)			○	○	OZB-A5703	
	1,0× (with micrometer) only in combination with OZB-A5703			○	○	OZB-A5704	
	for SLR cameras (Nikon)			○	○	OZB-A5706	
	for SLR cameras (Olympus)			○	○	OZB-A5707	
	for SLR cameras (Canon)			○	○	OZB-A5708	
Darkfield unit	Darkfield unit	○	○	○	○	OZB-A4601	
Object clamp	Object clamp	○	○	○	○	OBZ-A6205	
Stand	Pillar style, without illumination	✓		✓			
	Pillar style, with 3 W LED illumination (transmitted + incident)		✓		✓		
	Please find more stands in the catalogue on page 78 and on our website <a href="http://www.kern-sohn.com">www.kern-sohn.com</a>						
Stage plate	Frosted glass/∅ 94,5 mm		✓		✓	OZB-A5192	
	Black-white/∅ 94,5 mm	✓	✓	✓	✓	OZB-A5191	
	Glass/∅ 94,5 mm		○		○	OZB-A5190	
Mechanical stage (Pre-assembling on request)	Stage size W×D 188×160 mm, Travel 76×65 mm, for incident and transmitted illumination	○	○	○	○	OZB-A5781	
	Stage size W×D 180×175 mm, Travel 100×86 mm, for incident illumination only	○	○	○	○	OZB-A5782	
External illumination	Please find the information about external illumination units in the catalogue from page 81 and on our website <a href="http://www.kern-sohn.com">www.kern-sohn.com</a>						

✓ = Included with delivery

○ = Option

 360°	<b>360° rotatable microscope head</b>	 FL-LED	<b>Fluorescence illumination for compound microscopes</b> With 3 W LED illumination and filter	 SD	<b>SD card</b> For data storage
 MONO	<b>Monocular Microscope</b> For the inspection with one eye	 PH	<b>Phase contrast unit</b> For a higher contrast	 SOFTWARE	<b>PC software</b> To transfer the measurements from the device to a PC.
 BINO	<b>Binocular Microscope</b> For the inspection with both eyes	 DF	<b>Darkfield condenser/unit</b> For a higher contrast due to indirect illumination	 ATC	<b>Automatic temperature compensation</b> For measurements between 10 °C and 30 °C
 TRINO	<b>Trinocular Microscope</b> For the inspection with both eyes and the additional option for the connection of a camera	 POLAR	<b>Polarising unit</b> To polarise the light	 IP	<b>Protection against dust and water splashes IPxx</b> The type of protection is shown by the pictogram.
 ABBE	<b>Abbe Condenser</b> With high numerical aperture for the concentration and the focusing of light	 INFINITY	<b>Infinity system</b> Infinity corrected optical system	 BATT	<b>Battery operation</b> Ready for battery operation. The battery type is specified for each device.
 HAL	<b>Halogen illumination</b> For pictures bright and rich in contrast	 ZOOM	<b>Zoom magnification</b> For stereomicroscopes	 RECHARGE	<b>Battery operation rechargeable</b> Prepared for a rechargeable battery operation
 LED	<b>LED illumination</b> Cold, energy saving and especially long-life illumination	 PARALLEL	<b>Parallel optical system</b> For stereomicroscopes, enables fatigue-proof working	 230 V	<b>Mains adapter</b> 230V/50Hz in standard version for EU. On request GB, AUS or USA version.
 IL	<b>Incident illumination</b> For non-transparent objects	 SCALE	<b>Integrated scale</b> In the eyepiece	 230 V	<b>Power supply</b> Integrated in balance. 230V/50Hz standard EU. More standards e.g. GB, AUS or USA on request.
 TL	<b>Transmitting illumination</b> For transparent objects	 USB 2.0	<b>Integrated USB 2.0 digital camera</b> For direct transmitting of the picture to a PC	 1 DAY	<b>Package shipment</b> The time required to manufacture the product internally is shown in days in the pictogram.
 FL	<b>Fluorescence illumination</b> For stereomicroscopes	 USB 3.0	<b>Integrated USB 3.0 digital camera</b> For direct transmitting of the picture to a PC	 3 YEARS WARRANTY	<b>Warranty</b> The warranty period is shown in the pictogram.
 FL-HBO	<b>Fluorescence illumination for compound microscopes</b> With 100 W mercury lamp and filter	 HDMI	<b>HDMI digital camera</b> For direct transmitting of the picture to a display device		

## Abbreviations

<b>C-Mount</b>	Adapter for the connection of a camera to a trinocular microscope	<b>LWD</b>	Long Working Distance	<b>SWF</b>	Super Wide Field (Field number at least Ø 23 mm for 10x eyepiece)
<b>FPS</b>	Frames per second	<b>N.A.</b>	Numerical Aperture	<b>W.D.</b>	Working Distance
<b>H(S)WF</b>	High (Super) Wide Field (Eyepiece with high eye point for wearers of glasses)	<b>SLR Kamera</b>	Single-Lens Reflex camera	<b>WF</b>	Wide Field (Field number up to Ø 22 mm for 10x eyepiece)

## Your KERN specialist dealer: